

QA3

Patient Trolley System

Designed following comprehensive research into patient trolley function and ergonomics



Introducing the **NEW**
QA3 **DRIVE**

QA3 Patient Trolley System
Introduction
Overview
QA3 Patient Trolley Features
QA3 Patient Trolley Build Options
QA3 Drive Patient Trolley Features
QA3 Drive Patient Trolley Build Options
QA3 Drive Powered KneeFlex Patient Trolley Features
QA3 Drive Powered KneeFlex Patient Trolley Build Options
Optional Extras For Patient Trolleys
QA3 Emergency Trolley Features
QA3 Emergency Trolley Build Options
QA3 Drive Emergency Trolley Features
QA3 Drive Emergency Trolley Build Options
QA3 Drive Powered Emergency Trolley Features
QA3 Drive Powered Emergency Trolley Build Options
Optional Extras For Emergency Trolleys
In The Treatment Room...
Summary
Contacts



NEW

QA3 Patient Trolley System

Whether for Theatre, Accident & Emergency or any other hospital department or clinical area, the QA3 Variable Height Patient Trolley System is a compact, light weight and manoeuvrable patient trolley designed to combine maximum patient comfort and maximum ease of use for practitioners.

New **QA3 // // // // DRIVE** options include an inbuilt motorised drive, powered functions and KneeFlex leg positioning – read on for more details.

QA3
Patient
Trolley



QA3 Patient Trolley

The QA3 version 3.0 Variable Height Patient Trolley: ideal for patient transport, examination and intubation in theatre and clinical areas. Designed for maximum patient comfort, the trolley is supremely manoeuvrable, compact, light weight and robust.



3 YEAR WARRANTY
Terms and Conditions apply

QA3
Emergency
Trolley



QA3 Emergency Trolley

A specially designed version of the QA3 for Accident & Emergency departments is available. This includes an X-ray translucent platform, a medirail and the K8 pressure care mattress. Optional features include a monitor shelf and a Venturi suction system.

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley
Features

QA3 Patient Trolley
Build Options

QA3 Drive Patient Trolley
Features

QA3 Drive Patient Trolley
Build Options

QA3 Drive Powered KneeFlex
Patient Trolley Features

QA3 Drive Powered KneeFlex
Patient Trolley Build Options

Optional Extras
For Patient Trolleys

QA3 Emergency Trolley
Features

QA3 Emergency Trolley
Build Options

QA3 Drive
Emergency Trolley Features

QA3 Drive
Emergency Trolley Build Options

QA3 Drive Powered
Emergency Trolley Features

QA3 Drive Powered
Emergency Trolley Build Options

Optional Extras
For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



AneticAid

Innovative technology – practically applied

Anetic Aid is a world leading designer and manufacturer of medical devices used extensively for patient transport, treatment and recovery.

Renowned for innovation and high quality we work closely with clinicians to continuously develop and improve our products to help medical teams provide the best experience possible for their patients.

Our offer includes a broad range of equipment including surgical and patient transport trolley systems, tourniquets, theatre furniture, operating table accessories and surgical instruments.



QA3
Patient Trolley

QA3//////DRIVE
Patient Trolley

QA3//////DRIVE
Powered ⚡ KneeFlex Patient Trolley

QA3
Emergency Trolley

QA3//////DRIVE
Emergency Trolley

QA3//////DRIVE
Powered ⚡ Emergency Trolley

With its commitment to using the latest technology and materials for maximum functionality, the company is proud to be launching the latest evolution of its QA3 Patient Trolley System.

New options include an inbuilt motorised drive and electronic positioning – including KneeFlex – and all models are ideal for patient transport, examination and intubation. There are also specialist versions available for radiography in theatre and emergency areas. All are designed for maximum ease of use combined with superior patient comfort.

The QA3's compact construction is light weight yet robust, and its advanced steering system makes it supremely manoeuvrable

and easily guided by one person: the new option for an inbuilt motorised drive operated by touch sensitive pushing handles is completely unique.

With the backrest raised, it is little more than 1.8m long and it descends lower and is narrower and shorter than nearly all other patient trolleys on the market. It descends to an exceptionally low height, while unique foldaway side rails store away under the mattress footprint to give a virtually zero transfer gap.

In non-powered models, the backrest is gas assisted for easy lifting and in all versions, the Trendelenburg (and reverse Trendelenburg) tilt can be operated instantly at any height from the head end.

QA3 Patient Trolley System
Introduction
Overview
QA3 Patient Trolley Features
QA3 Patient Trolley Build Options
QA3 Drive Patient Trolley Features
QA3 Drive Patient Trolley Build Options
QA3 Drive Powered KneeFlex Patient Trolley Features
QA3 Drive Powered KneeFlex Patient Trolley Build Options
Optional Extras For Patient Trolleys
QA3 Emergency Trolley Features
QA3 Emergency Trolley Build Options
QA3 Drive Emergency Trolley Features
QA3 Drive Emergency Trolley Build Options
QA3 Drive Powered Emergency Trolley Features
QA3 Drive Powered Emergency Trolley Build Options
Optional Extras For Emergency Trolleys
In The Treatment Room...
Summary
Contacts

QA3

Patient Trolley

Code: 21110

Trendelenburg tilt ① extremely smooth and controlled mechanism, allows rapid Trendelenburg to be achieved even when trolley is at its lowest height



Pushing handles ⑨ four ergonomically positioned pushing handles give supreme manoeuvrability either with the patient lying flat or sitting up



FEATURES

The QA3 Patient Trolley (Code: 21110): ideal for patient transport, examination and intubation in theatre and clinical areas.

Fixed transfusion pole ② innovative quick release safety locking mechanism, integrates into either head or foot end of side rails for easy access and storage

Vertically operated side rails ③ store away under the mattress giving a virtually zero transfer gap. Spring counterbalancing allows effortless ease of operation

Hydraulic lift mechanism ④ exceptional low height capability, adjusted by means of independently operating pedals located on both sides of the trolley

Castors ⑤ polyurethane tyres offer low rolling resistance for high manoeuvrability, 150mm castors fitted as standard, 200mm optional

5th wheel steering ⑥ unique mechanism which ensures trolley is highly manoeuvrable and easily guided by one person

Gas assisted backlift ⑦ unique design, precision controlled, easy lift system for enhanced user and patient comfort

Cylinder trough ⑧ side loading for ease of access, suitable for a wide range of cylinder options including F, E or D



QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



AneticAid

QA3

Patient Trolley

Code: 21110

BUILD OPTIONS

The QA3 Patient Trolley (Code: 21110) has a number of build options for specific clinical applications. Please specify at the time of ordering.



Factory Fitted Build Options The following Build Options must be specified at the time of ordering.

Code	Name
21130/3	Medirail – T-bar & V-bracket – QA3 v3.0 for Serres System
21130/4	Medirail – T-bar & V-bracket – QA3 v3.0 for Carnet System
21149/2	V-bracket for VacSax Canister
21132	200mm Castors
21135	Reverse Trendelenburg
21138	Increased Maximum Height 100mm
21155	Fixed Pushing Handles Foot End
21159	Paper Roll Holder (up to 620mm wide) – Universal

Code	Name
21161	Additional Transfusion Pole – Fixed
21175	Department Identification
21191	Patient's Storage Container
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
53556A	Oxygen Cylinder Strap for CD & E Sized Cylinders
53556A	Oxygen Cylinder Strap for F, HD & ZX Sized Cylinders
22186	K8 Mattress – with Zip Cover – Black
22186-C1	K8 Mattress – with Zip Cover – Internal Foam Protector Fitted



QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



QA3 DRIVE

Patient Trolley

Code: 21112

Trendelenburg tilt ① extremely smooth and controlled mechanism, allows rapid Trendelenburg to be achieved even when trolley is at its lowest height



NEW Pushing handles ⑪ totally innovative touch sensitive design to activate inbuilt motorised drive and move trolley forward or back



FEATURES

The QA3 DRIVE Patient Trolley (Code: 21112) has a new inbuilt motorised drive and is supremely manoeuvrable, compact, light weight and robust.

NEW Inbuilt motorised drive ② works with 5th wheel steering to offer superb manoeuvrability and controlled, effortless travel up to gradients of 10°

Fixed transfusion pole ③ innovative quick release safety locking mechanism, integrates into either head or foot end of side rail for easy access and storage

Vertically operated side rails ④ store away under the mattress giving a virtually zero transfer gap. Spring counterbalancing allows effortless ease of operation

Castors ⑤ polyurethane tyres offer low rolling resistance for high manoeuvrability, 150mm castors fitted as standard, 200mm optional

Gas assisted backlift ⑥ unique design, precision controlled, easy lift system for enhanced user and patient comfort

Cylinder trough ⑦ side loading for ease of access, suitable for a wide range of cylinder options including F, E or D

Infection control ⑧ smooth moulded surfaces and acrylic capped vac formings enclose the trolley mechanics to minimise dirt traps. Anti-bacterial properties and innovative design ensure the product is easy to clean

K8 pressure care mattress ⑨ manufactured from temperature-sensitive gel feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Developed with the flexible properties of lycra, this outer material is anti-microbacterial, water resistant, vapour permeable and latex free

Patient comfort ⑩ exceptional low height capability and store away side rails allow easy mount/dismount particularly for elderly patients. Wider patient platform and pressure care mattress ensure maximum patient comfort. Improved ground clearance allows better access for patient hoists



MADE IN THE UK QUALITY ASSURED GLOBAL DISTRIBUTION



QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



AneticAid

QA3 DRIVE

Patient Trolley

Code: 21112

BUILD OPTIONS


The QA3 DRIVE Patient Trolley (Code: 21112) has a number of build options for specific clinical applications. Please specify at the time of ordering.



21149/2



21159

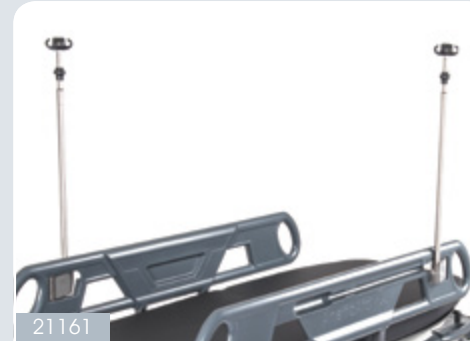

21161



21191



21393



22186

Factory Fitted Build Options The following Build Options must be specified at the time of ordering.

Code	Name
21130/3	Medirail – T-bar & V-bracket – QA3 v3.0 for Serres System
21130/4	Medirail – T-bar & V-bracket – QA3 v3.0 for Carnet System
21149/2	V-bracket for VacSax Canister
21159	Paper Roll Holder (up to 620mm wide) – Universal
21161	Additional Transfusion Pole – Fixed
21175	Department Identification

Code	Name
21191	Patient's Storage Container
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
53556A	Oxygen Cylinder Strap for CD & E Sized Cylinders
53556A	Oxygen Cylinder Strap for F, HD & ZX Sized Cylinders
22186	K8 Mattress – with Zip Cover – Black
22186-C1	K8 Mattress – with Zip Cover – Internal Foam Protector Fitted



QUALITY UK PRODUCTION

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts

QA3 DRIVE

Powered  KneeFlex Patient Trolley

Code: 21114

NEW Positive feel push button controls **1** for electrically powered functions



NEW CPR button **2** automatically positions trolley flat for CPR



NEW KneeFlex **12** smooth electrically powered KneeFlex position for additional patient comfort

FEATURES

The QA3 DRIVE Powered KneeFlex Patient Trolley (Code: 21112) has all the benefits of the QA3 Drive – plus electrically powered functions including a 'KneeFlex' facility for additional patient comfort

NEW Inbuilt motorised drive **3** works with 5th wheel steering to offer superb manoeuvrability and controlled, effortless travel up to gradients of 10°

NEW Pushing handles **4** totally innovative touch sensitive design to activate inbuilt motorised drive and move trolley forward or back

NEW Electrically powered backlift **5** unique design, precision controlled, for enhanced user and patient comfort

NEW Electrically powered raise and lower mechanism **6** exceptional low height capability

Trendelenburg tilt **7** extremely smooth and controlled mechanism, allows rapid Trendelenburg to be achieved even when trolley is at its lowest height

Vertically operated side rails **8** store away under the mattress giving a virtually zero transfer gap. Spring counterbalancing allows effortless ease of operation

Fixed transfusion pole **9** innovative quick release safety locking mechanism, integrates into either head or foot end of side rails for easy access and storage

Cylinder trough **10** side loading for ease of access, suitable for a wide range of cylinder options including F, E or D

Castors **11** polyurethane tyres offer low rolling resistance for high manoeuvrability, 150mm castors fitted as standard, 200mm optional



MADE IN THE UK



QUALITY ASSURED



GLOBAL DISTRIBUTION



QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



AneticAid

QA3 DRIVE

Powered  KneeFlex Patient Trolley

Code: 21114

BUILD OPTIONS

The QA3 DRIVE Powered KneeFlex Patient Trolley (Code: 21114) has a number of build options for specific clinical applications. Please specify at the time of ordering.



21149/2



21161



21175



21191



21393



22186

Factory Fitted Build Options The following Build Options must be specified at the time of ordering.

Code	Name
21130/3	Medirail – T-bar & V-bracket – QA3 v3.0 for Serres System
21130/4	Medirail – T-bar & V-bracket – QA3 v3.0 for Carnet System
21149/2	V-bracket for VacSax Canister
21161	Additional Transfusion Pole – Fixed
21175	Department Identification
21191	Patient's Storage Container

Code	Name
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
53556A	Oxygen Cylinder Strap for CD & E Sized Cylinders
53556A	Oxygen Cylinder Strap for F, HD & ZX Sized Cylinders
22186	K8 Mattress – with Zip Cover – Black
22186-C1	K8 Mattress – with Zip Cover – Internal Foam Protector Fitted



QUALITY UK PRODUCTION

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts

Optional
Extras For:

QA3
Patient Trolley
Code: 21110

QA3 // // // DRIVE
Patient Trolley
Code: 21112

QA3 // // // DRIVE
Powered ⚡ KneeFlex Patient Trolley
Code: 21114



21145



21152-A1



21393



21173



21188



21395

The following Optional Extras can be purchased and used with the trolley at any time.

Code	Name
21143	Suction – Injector & Probe for CD and HX Cylinder – Complete
21145	Suction – Emergency Venturi System – Complete
21147	1 Litre VacSax Standard Bore Liners (Box of 25)
21152	Monitor Shelf – with Removable Refreshment Tray
21152-B1	Monitor Shelf – Transfusion Poles (Factory Fitted Option)
21152-A1	Refreshment Tray Only
21157	Lateral Arm Positioner with 40mm Pad

Code	Name
21165	Endoscopy Height Extension Demi-frame
21172	Side Rail – Grey Covers – Pair
21173	Side Rail – Child Print Covers – Pair
21188	Theatre Pillow with Mattress Retention Strap
21390	Push / Pull Bar
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
21395	Foot-end Extension c/w 30mm Pad

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley
Features

QA3 Patient Trolley
Build Options

QA3 Drive Patient Trolley
Features

QA3 Drive Patient Trolley
Build Options

QA3 Drive Powered KneeFlex
Patient Trolley Features

QA3 Drive Powered KneeFlex
Patient Trolley Build Options

Optional Extras
For Patient Trolleys

QA3 Emergency Trolley
Features

QA3 Emergency Trolley
Build Options

QA3 Drive
Emergency Trolley Features

QA3 Drive
Emergency Trolley Build Options

QA3 Drive Powered
Emergency Trolley Features

QA3 Drive Powered
Emergency Trolley Build Options

Optional Extras
For Emergency Trolleys

In The Treatment Room...


Summary

Contacts

QA3

Emergency Trolley

Code: 21120

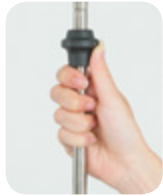

2 fixed transfusion poles ① fitted as standard to emergency trolley



Patient comfort ⑥ exceptional low height capability and store away side rails allow easy mount/dismount particularly for elderly patients. Wider patient platform and pressure care mattress ensure maximum patient comfort. Improved ground clearance allows better access for patient hoists



FEATURES

The QA3 Drive Emergency Trolley (Code: 21120) is specially designed for Accident & Emergency departments. It includes an X-ray translucent platform, a medirail and the K8 pressure care mattress.

X-ray translucent platform ② full length radio translucent platform suitable for digital or traditional plates

Infection control ③ smooth moulded surfaces and acrylic capped vac formings enclose the trolley mechanics to minimise dirt traps. Anti-bacterial properties and innovative design ensure the product is easy to clean

Advanced braking ④ exceptionally stable platform, advanced single bar braking system that allows all four castors to be braked simultaneously from either side of the trolley

K8 pressure care mattress ⑤ manufactured from temperature-sensitive gel feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Developed with the flexible properties of lycra, this outer material is anti-microbacterial, water resistant, vapour permeable and latex free

K8

TECHNOLOGY



MADE IN THE UK QUALITY ASSURED GLOBAL DISTRIBUTION



QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



AneticAid

QA3

Emergency Trolley

Code: 21120

BUILD OPTIONS

The QA3 Emergency Trolley (Code: 21120) has a number of build options for specific clinical applications. Please specify at the time of ordering.



21130/3



21132



21135



21175



21191



22186

Factory Fitted Build Options The following Build Options must be specified at the time of ordering.

Code	Name
21130/3	Medirail – T-bar & V-bracket – QA3 v3.0 for Serres System
21130/4	Medirail – T-bar & V-bracket – QA3 v3.0 for Carnet System
21132	200mm Castors
21135	Reverse Trendelenburg
21138	Increased Maximum Height 100mm
21155	Fixed Pushing Handles Foot End
21159	Paper Roll Holder (up to 620mm wide) – Universal
21175	Department Identification
21191	Patient's Storage Container

Code	Name
21196	X-ray Tray Guides
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
53556A	Oxygen Cylinder Strap for CD & E Sized Cylinders
53556A	Oxygen Cylinder Strap for F, HD & ZX Sized Cylinders
22134	5th Wheel Steer – Extended Activation Pedal
TBC	5th Wheel Steer – Foot-end Mounted
22186	K8 Mattress – with Zip Cover – Black
22186-C1	K8 Mattress – with Zip Cover – Internal Foam Protector Fitted



QUALITY UK PRODUCTION

- QA3 Patient Trolley System
- Introduction
- Overview
- QA3 Patient Trolley Features
- QA3 Patient Trolley Build Options
- QA3 Drive Patient Trolley Features
- QA3 Drive Patient Trolley Build Options
- QA3 Drive Powered KneeFlex Patient Trolley Features
- QA3 Drive Powered KneeFlex Patient Trolley Build Options
- Optional Extras For Patient Trolleys
- QA3 Emergency Trolley Features
- QA3 Emergency Trolley Build Options**
- QA3 Drive Emergency Trolley Features
- QA3 Drive Emergency Trolley Build Options
- QA3 Drive Powered Emergency Trolley Features
- QA3 Drive Powered Emergency Trolley Build Options
- Optional Extras For Emergency Trolleys
- In The Treatment Room...
- Summary
- Contacts

QA3 DRIVE

Emergency Trolley

Code: 21122



2 fixed transfusion poles ① fitted as standard to emergency trolley

Patient comfort ⑪ exceptional low height capability and store away side rails allow easy mount/dismount particularly for elderly patients. Wider patient platform and pressure care mattress ensure maximum patient comfort. Improved ground clearance allows better access for patient hoists



FEATURES

The QA3 DRIVE Emergency Trolley (Code: 21122) has an inbuilt motorised drive and also offers an X-ray translucent platform, a medirail and the K8 pressure care mattress.

NEW Inbuilt motorised drive ② works with 5th wheel steering to offer superb manoeuvrability and controlled, effortless travel up to gradients of 10°

NEW Pushing handles ③ totally innovative touch sensitive design to activate inbuilt motorised drive and move trolley forward or back

Trendelenburg tilt ④ extremely smooth and controlled mechanism, allows rapid Trendelenburg to be achieved even when trolley is at its lowest height

Vertically operated side rails ⑤ store away under the mattress giving a virtually zero transfer gap. Spring counterbalancing allows effortless ease of operation

Cylinder trough ⑥ side loading for ease of access, suitable for a wide range of cylinder options including F, E or D

Castors ⑦ polyurethane tyres offer low rolling resistance for high manoeuvrability, 150mm castors fitted as standard, 200mm optional

Infection control ⑧ smooth moulded surfaces and acrylic capped vac formings enclose the trolley mechanics to minimise dirt traps. Anti-bacterial properties and innovative design ensure the product is easy to clean

K8 pressure care mattress ⑨ manufactured from temperature-sensitive gel feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Developed with the flexible properties of lycra, this outer material is anti-microbial, water resistant, vapour permeable and latex free

X-ray translucent platform ⑩ full length radio translucent platform suitable for digital or traditional plates



QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts

QA3 DRIVE

Emergency Trolley

Code: 21122

BUILD OPTIONS

The QA3 DRIVE Emergency Trolley (Code: 21120) has a number of build options for specific clinical applications. Please specify at the time of ordering.



21130/3



21175



21191



21196



21393



22186

Factory Fitted Build Options The following Build Options must be specified at the time of ordering.

Code	Name
21130/3	Medirail – T-bar & V-bracket – QA3 v3.0 for Serres System
21130/4	Medirail – T-bar & V-bracket – QA3 v3.0 for Carnet System
21175	Department Identification
21191	Patient's Storage Container
21196	X-ray Tray Guides

Code	Name
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
53556A	Oxygen Cylinder Strap for CD & E Sized Cylinders
53556A	Oxygen Cylinder Strap for F, HD & ZX Sized Cylinders
22186	K8 Mattress – with Zip Cover – Black
22186-C1	K8 Mattress – with Zip Cover – Internal Foam Protector Fitted



QUALITY UK PRODUCTION

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts

QA3 DRIVE

Powered ⚡ Emergency Trolley

Code: 21124

NEW Positive feel push button controls ① for electrically powered functions



Patient comfort ⑫ exceptional low height capability and store away side rails allow easy mount/dismount particularly for elderly patients. Wider patient platform and pressure care mattress ensure maximum patient comfort. Improved ground clearance allows better access for patient hoists



FEATURES

The QA3 DRIVE Powered Emergency Trolley (Code: 21124) offers inbuilt motorised drive and powered patient positioning: raising and lowering the whole trolley platform and the backrest.

NEW Pushing handles ② totally innovative touch sensitive design to activate inbuilt motorised drive and move trolley forward or back

NEW Inbuilt motorised drive ③ works with 5th wheel steering to offer superb manoeuvrability and controlled, effortless travel up to gradients of 10°

Trendelenburg tilt ④ extremely smooth and controlled mechanism, allows rapid Trendelenburg to be achieved even when trolley is at its lowest height

Vertically operated side rails ⑤ store away under the mattress giving a virtually zero transfer gap. Spring counterbalancing allows effortless ease of operation

2 fixed transfusion poles ⑥ fitted as standard to emergency trolley

Cylinder trough ⑦ side loading for ease of access, suitable for a wide range of cylinder options including F, E or D

Castors ⑧ polyurethane tyres offer low rolling resistance for high manoeuvrability, 150mm castors fitted as standard, 200mm optional

Infection control ⑨ smooth moulded surfaces and acrylic capped vac formings enclose the trolley mechanics to minimise dirt traps. Anti-bacterial properties and innovative design ensure the product is easy to clean

K8 pressure care mattress ⑩ manufactured from temperature-sensitive gel feel foam that moulds to the patient's body shape, relieving pressure by distributing the load as evenly as possible. Developed with the flexible properties of lycra, this outer material is anti-microbacterial, water resistant, vapour permeable and latex free

X-ray translucent platform ⑪ full length radio translucent platform suitable for digital or traditional plates



- QA3 Patient Trolley System
- Introduction
- Overview
- QA3 Patient Trolley Features
- QA3 Patient Trolley Build Options
- QA3 Drive Patient Trolley Features
- QA3 Drive Patient Trolley Build Options
- QA3 Drive Powered KneeFlex Patient Trolley Features
- QA3 Drive Powered KneeFlex Patient Trolley Build Options
- Optional Extras For Patient Trolleys
- QA3 Emergency Trolley Features
- QA3 Emergency Trolley Build Options
- QA3 Drive Emergency Trolley Features
- QA3 Drive Emergency Trolley Build Options
- QA3 Drive Powered Emergency Trolley Features**
- QA3 Drive Powered Emergency Trolley Build Options
- Optional Extras For Emergency Trolleys
- In The Treatment Room...
- Summary
- Contacts



QA3 DRIVE

Powered  Emergency Trolley

Code: 21124

BUILD OPTIONS

The QA3 DRIVE Powered Emergency Trolley (Code: 21124) has a number of build options for specific clinical applications. Please specify at the time of ordering.



21130/3



21175



21191



21196



21393



22186

Factory Fitted Build Options The following Build Options must be specified at the time of ordering.

Code	Name
21130/3	Medirail – T-bar & V-bracket – QA3 v3.0 for Serres System
21130/4	Medirail – T-bar & V-bracket – QA3 v3.0 for Carnet System
21175	Department Identification
21191	Patient's Storage Container
21196	X-ray Tray Guides

Code	Name
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
53556A	Oxygen Cylinder Strap for CD & E Sized Cylinders
53556A	Oxygen Cylinder Strap for F, HD & ZX Sized Cylinders
22186	K8 Mattress – with Zip Cover – Black
22186-C1	K8 Mattress – with Zip Cover – Internal Foam Protector Fitted



QUALITY UK PRODUCTION

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts

Optional
Extras For:

QA3
Emergency Trolley
Code: 21120

QA3 // DRIVE
Emergency Trolley
Code: 21122

QA3 // DRIVE
Powered ⚡ Emergency Trolley
Code: 21114



21145



21152



21152-A1



21158



21189



21192

The following Optional Extras can be purchased and used with the trolley at any time.

Code	Name
21143	Suction – Injector & Probe for CD and HX Cylinder – Complete
21145	Suction – Emergency Venturi System – Complete
21147	1 Litre VacSax Standard Bore Liners (Box of 25)
21152	Monitor Shelf – with Removable Refreshment Tray
21152-B1	Monitor Shelf – Transfusion Poles *Factory Fitted Option
21152-A1	Refreshment Tray Only
21153	Foot-end Extension with 100mm Pad
21158	Lateral Arm Positioner with 100mm Pad

Code	Name
21172	Side Rail – Grey Covers – Pair
21173	Side Rail – Child Print Covers – Pair
21188	Theatre Pillow with Mattress Retention Strap
21189	Foot-end Mattress Protector
21192	Lateral Cassette Holder
21195	X-ray Tray
21390	Push / Pull Bar

- QA3 Patient Trolley System
- Introduction
- Overview
- QA3 Patient Trolley Features
- QA3 Patient Trolley Build Options
- QA3 Drive Patient Trolley Features
- QA3 Drive Patient Trolley Build Options
- QA3 Drive Powered KneeFlex Patient Trolley Features
- QA3 Drive Powered KneeFlex Patient Trolley Build Options
- Optional Extras For Patient Trolleys
- QA3 Emergency Trolley Features
- QA3 Emergency Trolley Build Options
- QA3 Drive Emergency Trolley Features
- QA3 Drive Emergency Trolley Build Options
- QA3 Drive Powered Emergency Trolley Features
- QA3 Drive Powered Emergency Trolley Build Options
- Optional Extras For Emergency Trolleys
- In The Treatment Room...
- Summary
- Contacts

In the Treatment Room...

All versions of the QA3 Patient Trolley are excellent options in treatment rooms such as Endoscopy Suites – and we recommend a number of specifications to tailor them for use in these areas.

QA3
Patient Trolley

Code: 21110

QA3//////**DRIVE**
Patient Trolley

Code: 21112

QA3//////**DRIVE**
Powered ⚡ KneeFlex Patient Trolley

Code: 21114

Factory Fitted Build Options

The following Build Options must be specified at the time of ordering.



QUALITY UK PRODUCTION




21135



21138



21175

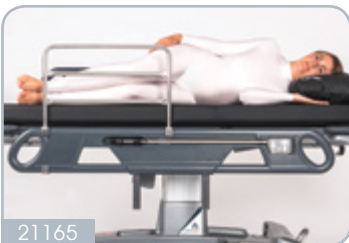

21165



21152-A1



21173

Additional Factory Fitted Build Options

Code	Name
21135	Reverse Trendelenburg
21138	Increased Maximum Height 100mm
21161	Additional Transfusion Pole – Fixed
21149/2	V-bracket for VacSax Canister
21132	200mm Castors
21175	Department Identification
21393	Oxygen Cylinder Bracket for CD Systems with Backing Plate
53556A	Oxygen Cylinder Strap for CD & E Sized Cylinders
53556A	Oxygen Cylinder Strap for F, HD & ZX Sized Cylinders
22186	K8 Mattress – with Zip Cover – Black
22186-C1	K8 Mattress – with Zip Cover – Internal Foam Protector Fitted

Optional Extras

Code	Name
21143	Suction – Injector & Probe for CD and HX Cylinder – Complete
21145	Suction – Emergency Venturi System – Complete
21147	1 Litre VacSax Standard Bore Liners (Box of 25)
21152-A1	Refreshment Tray Only
21158	Lateral Arm Positioner with 100mm Pad
21165	Endoscopy Height Extension Demi-frame
21172	Side Rail – Grey Covers – Pair
21173	Side Rail – Child Print Covers – Pair
21188	Theatre Pillow with Mattress Retention Strap
21390	Push / Pull Bar
21395	Foot-end Extension c/w 30mm Pad

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



AneticAid

Summary

ITEM	Cat No.
QA3 Patient Trolley v3.0	21110
QA3 DRIVE Patient Trolley	21112
QA3 DRIVE Powered KneeFlex Patient Trolley	21114
QA3 Emergency Trolley v3.0 with Radiolucent Platform	21120
QA3 DRIVE Emergency Trolley with Radiolucent Platform	21122
QA3 DRIVE Powered Emergency Trolley with Radiolucent Platform	21124

DIMENSIONS

Overall length: 2130mm

Overall width: 735mm

Height range (Patient Trolley): 475mm to 775mm (excluding mattress)

Height range (Emergency Trolley): 520mm to 820mm (excluding mattress)

Trolley weight: 128kg (21110) 148kg (21120)

Maximum patient weight: 320kg (Please refer to user manual)

QA3 Patient Trolley System
Introduction
Overview
QA3 Patient Trolley Features
QA3 Patient Trolley Build Options
QA3 Drive Patient Trolley Features
QA3 Drive Patient Trolley Build Options
QA3 Drive Powered KneeFlex Patient Trolley Features
QA3 Drive Powered KneeFlex Patient Trolley Build Options
Optional Extras For Patient Trolleys
QA3 Emergency Trolley Features
QA3 Emergency Trolley Build Options
QA3 Drive Emergency Trolley Features
QA3 Drive Emergency Trolley Build Options
QA3 Drive Powered Emergency Trolley Features
QA3 Drive Powered Emergency Trolley Build Options
Optional Extras For Emergency Trolleys
In The Treatment Room...
Summary
Contacts

QA3

Patient Trolley System



For more information on any of our products or service contracts, call:

T +44 (0)1943 878647

sales@aneticaid.com

aneticaid.com

SPECIALISING IN:

Patient & Surgery Trolley Systems

Operation Table Accessories

Stainless Steel Theatre Furniture

Surgical Instruments

Electro-Surgical Accessories

Tourniquet Systems

Fibre Optic Instruments

Service & Maintenance

Service and Maintenance

All Anetic Aid products are fully supported and it is our policy to work closely with our distributors to ensure each piece of equipment is maintained to the highest standard.

Distributor:

For details of your local distributor of Anetic Aid products, call **T** +44 (0)1943 878647



MADE IN THE UK



QUALITY ASSURED



GLOBAL DISTRIBUTION

QA3 Patient Trolley System

Introduction

Overview

QA3 Patient Trolley Features

QA3 Patient Trolley Build Options

QA3 Drive Patient Trolley Features

QA3 Drive Patient Trolley Build Options

QA3 Drive Powered KneeFlex Patient Trolley Features

QA3 Drive Powered KneeFlex Patient Trolley Build Options

Optional Extras For Patient Trolleys

QA3 Emergency Trolley Features

QA3 Emergency Trolley Build Options

QA3 Drive Emergency Trolley Features

QA3 Drive Emergency Trolley Build Options

QA3 Drive Powered Emergency Trolley Features

QA3 Drive Powered Emergency Trolley Build Options

Optional Extras For Emergency Trolleys

In The Treatment Room...

Summary

Contacts



Anetic Aid Ltd. Sapper Jordan Rossi Business Park, Otley Road, Baildon BD17 7AX, UK